


ORIGIO® Handling™ Medium



An advanced medium providing a stable pH environment for handling and micromanipulation of oocytes outside the CO₂ incubator

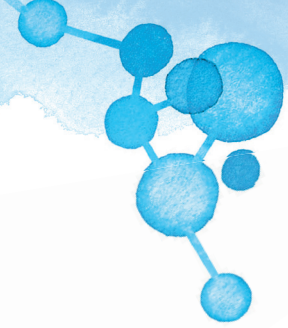


REF 83110060D
ORIGIO® Handling™
with phenol red
Exp. date: 2020-02-02
Lot: 12345678
Store at 2-8°C. Protect from light
Use within 7 days after opening.
Discard excess (unused) media
following warming. Aseptically filtered
Info: See package insert. Rx only
Date opened: _____
origio 60 mL
ORIGIO as DK-2780 Målev, Denmark

Indication for use: For in vitro procedures involving handling and micromanipulation of gametes and embryos outside the CO₂ incubator
Indication includes oocyte retrieval including follicle flushing, gamete (oocytes& sperm), embryo washing, micromanipulation procedures including ICSI, assisted hatching and trophoctoderm biopsy
Endotoxin tested 50.1 EU/ml
1 cell MEA 2: 80% blastocysts at 96h
HSSA: Sperm Motility 2: 86% of control after 6h
Contains Gentamicin 10 µg/ml and HSA (US licensed source)



CooperSurgical®
Fertility Solutions

ORIGIO[®] Handling[™] Medium

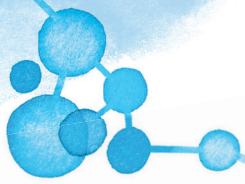

ORIGIO Handling is a modern, ready-to-use handling medium with a dual buffer (HEPES/MOPS) that meets the need for stable pH level outside of the incubator.

HEPES and MOPS, bicarbonate buffered

Optimized external pH buffer capacity

Essential amino acids

Support optimal oocyte viability/development

Vitamins

Act as antioxidants to reduce cellular stress

Physiological composition in line with ORIGIO Sequential Fert[™]

Complementary to the ORIGIO Sequential Series

Product Specifications

- HEPES / MOPS
- HSA: 5 mg/mL
- Available with and without phenol red
- Shelf life: 26 weeks from date of manufacture
- Shelf life after opening: 7 days
- Ready-to-use

Components

- Physiological salts
- HEPES / MOPS
- Bicarbonate
- Energy substrates
- Amino acids
- Vitamins
- HSA
- Gentamicin
- Phenol red

Product Code	Product Name	Volume (mL)	Description
83100060	ORIGIO Handling	60	Without phenol red
83100125	ORIGIO Handling	125	Without phenol red
83110060	ORIGIO Handling	60	With phenol red
83110125	ORIGIO Handling	125	With phenol red

CE
2797

